

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

1634
17w

Applicant: Jiri Snaidr

Title: METHOD OF DETECTING MICROORGANISMS IN A SAMPLE

Docket No.: 235.017US1

Filed: November 7, 2001

Examiner: Sally A. Sakelaris



Serial No.: 10/008,523

Due Date: N/A

Group Art Unit: 1634

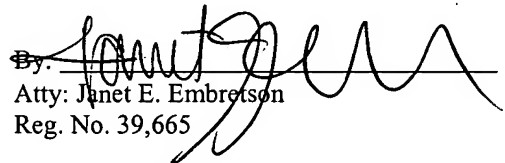
Commissioner for Patents
P.O. Box 1450
Alexandria, VA 22313-1450

We are transmitting herewith the following attached items (as indicated with an "X"):

- ☒ A return postcard.
- ☒ Declaration Under 37 CFR 1.132 (2 Pages).

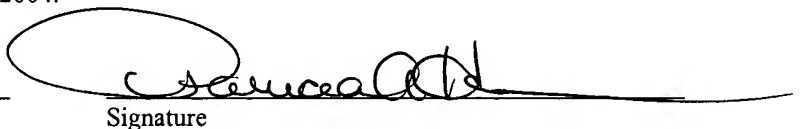
If not provided for in a separate paper filed herewith, Please consider this a PETITION FOR EXTENSION OF TIME for sufficient number of months to enter these papers and please charge any additional fees or credit overpayment to Deposit Account No. 19-0743.

SCHWEGMAN, LUNDBERG, WOESSNER & KLUTH, P.A.
Customer Number 21186

By: 
Atty: Janet E. Embretson
Reg. No. 39,665

CERTIFICATE UNDER 37 CFR 1.8: The undersigned hereby certifies that this correspondence is being deposited with the United States Postal Service with sufficient postage as first class mail, in an envelope addressed to: Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450, on this 3rd day of June, 2004.

PATRICIA A. HULTMAN
Name


Signature

SCHWEGMAN, LUNDBERG, WOESSNER & KLUTH, P.A.
(GENERAL)

**S/N 10/008,523****PATENT****IN THE UNITED STATES PATENT AND TRADEMARK OFFICE**

Applicant:	Jiri Snaidr	Examiner:	Sally A. Sakelaris
Serial No.:	10/008,523	Group Art Unit:	1634
Filed:	November 7, 2001	Docket No.:	235.017US1
Title:	METHOD OF DETECTING MICROORGANISMS IN A SAMPLE		

DECLARATION UNDER 37 C.F.R. § 1.132

Commissioner for Patents
P.O. Box 1450
Alexandria, VA 22313-1450

Sir:

I, Dr. Jiri Snaidr, declare as follows:

1. I am the inventor of the claims of the present application and make this Declaration in support of the patentability of the claims of this application as amended in the Amendment which accompanies this Declaration.
2. The claims of the above-identified application are directed to methods to detect microorganisms which employ a particular separation solution, e.g., one having water, 0.001-0.01 M Tris/HCl, pH 9.0 +/- 2.0, DMSO or 1X SSC.
3. As described at page 8, line 25-page 9, line 13 of the specification, experiments were performed with the following separation solutions at 80°C: 0.01 M Tris-HCl, pH 9.0, water, 1 X SSC, pH 10.0, DMSO, and formamide; and at 100°C: 0.01 M Tris-HCl, pH 9.0, water, and formamide. It is disclosed that all of the non-formamide separation solutions provide a better signal than a formamide separation solution (page 8, line 34-page 9, line 3 and page 9, lines 6-11).
4. In addition, Cy3-labeled oligonucleotides were diluted in water or formamide and incubated at 80°C for up to 15 minutes (Figure 1). The fluorescent signal obtained from oligonucleotides diluted in formamide was about 80% lower than the signal obtained from oligonucleotides diluted in water.



Declaration Under 37 CFR 1.132

Serial No.: 10/008,523

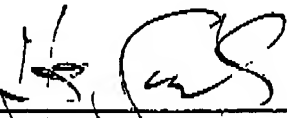
Filed: November 7, 2001

Title: METHOD OF DETECTING MICROORGANISMS IN A SAMPLE

Page 2 of 2

SLWK 235.017US1

5. Moreover, the signal from Cy3-labeled oligonucleotides diluted in 0.01 M Tris-HCl, pH 9, and incubated at 80°C was similar to the signal obtained with DMSO (Figure 2).
6. Thus, the use of separation solutions within the scope of the claims to denature target hybridized probes having a detectable signal yields a stronger signal than the use of a formamide separation solution to denature those probes.
7. I hereby declare that all statements made herein of our own knowledge are true and that all statements made on information and belief are believed to be true; and further that the statements were made with the knowledge that willful false statement and the like so made are punishable by fine or imprisonment, or both, under Section 1001 of Title 18 of the United States Code, and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

Dated: June 1, 2004By: 
Dr. Jiri Snajdr